


# Manual EZwash™

Low cost high performance

EZwash is the latest convenient attachment available for the Multi-tool handset, for the manual cleaning of microplates. The EZwash is compact, easy to use, cost effective and autoclavable. Comes complete with stand that suits 8 and 12 way manifolds with built-in height adjustment for different plate depths.



The EZwash uses gravitational fill via the angled stainless steel tubes and vacuums fluid from the other. However it's also possible to use with two vacuum systems to give powered fill and vacuum.

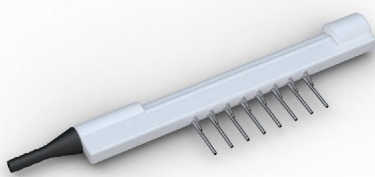

8way manifold



12way manifold

\* Zisys Ltd reserve the right to change the design without prior notice as part of the ongoing design process.



Precision Vacuum System



Multi-tool Aspiration Handset

- 8way EZwash complete kit
- 12way EZwash complete kit
- 8way EZwash
- 12way EZwash

**zisys**<sup>TM</sup>

39 Mackintosh Place  
South Newmoor Industrial Estate  
Irvine  
Ayrshire  
Scotland KA11 4JT  
Tel: +44 (0) 1294 219266  
[www.zisysautomation.com](http://www.zisysautomation.com)  
[allan@zisysautomation.com](mailto:allan@zisysautomation.com)

\* Zisys Ltd reserve the right to change the design without prior notice as part of the ongoing design process.